



Designation: F2004 – 16

Standard Test Method for Transformation Temperature of Nickel-Titanium Alloys by Thermal Analysis¹

This standard is issued under the fixed designation F2004; the number immediately following the designation indicates the year of original adoption or, in the case of revision, the year of last revision. A number in parentheses indicates the year of last reapproval. A superscript epsilon (ϵ) indicates an editorial change since the last revision or reapproval.

1. Scope

1.1 This test method defines procedures for determining the transformation temperatures of nickel-titanium shape memory alloys, produced in accordance with Specification F2063, by differential scanning calorimetry.

1.2 The values stated in SI units are to be regarded as standard. No other units of measurement are included in this standard.

1.3 *This standard does not purport to address all of the safety concerns, if any, associated with its use. It is the responsibility of the user of this standard to establish appropriate safety and health practices and to determine the applicability of regulatory limitations prior to use.*

2. Referenced Documents

2.1 ASTM Standards:²

- E177 Practice for Use of the Terms Precision and Bias in ASTM Test Methods
- E473 Terminology Relating to Thermal Analysis and Rheology
- E691 Practice for Conducting an Interlaboratory Study to Determine the Precision of a Test Method
- E967 Test Method for Temperature Calibration of Differential Scanning Calorimeters and Differential Thermal Analyzers
- E1142 Terminology Relating to Thermophysical Properties
- F2005 Terminology for Nickel-Titanium Shape Memory Alloys
- F2063 Specification for Wrought Nickel-Titanium Shape Memory Alloys for Medical Devices and Surgical Implants
- F2082 Test Method for Determination of Transformation

Temperature of Nickel-Titanium Shape Memory Alloys by Bend and Free Recovery

3. Terminology

3.1 Specific technical terms used in this test method are found in Terminologies E473, E1142, and F2005.

4. Summary of Test Method

4.1 This test method involves heating and cooling a test specimen at a controlled rate in a controlled environment through the temperature interval of the phase transformation. The difference in heat flow between the test material and a reference material due to energy changes is continuously monitored and recorded. Absorption of energy due to a phase transformation in the specimen results in an endothermic peak on heating. Release of energy due to a phase transformation in the specimen results in an exothermic peak on cooling.

5. Significance and Use

5.1 Differential scanning calorimetry provides a rapid method for determining the transformation temperature(s) of nickel-titanium shape memory alloys.

5.2 This test method uses small, stress-free, annealed samples to determine whether a sample of nickel-titanium alloy containing nominally 54.5 to 56.5 % nickel by weight is austenitic or martensitic at a particular temperature. Since chemical analysis of these alloys does not have sufficient precision to determine the transformation temperature by measuring the nickel-to-titanium ratio of the alloy, direct measurement of the transformation temperature of an annealed sample of known thermal history is recommended.

5.3 This test method is useful for quality control, specification acceptance, and research.

5.4 Transformation temperatures derived from differential scanning calorimetry (DSC) may not agree with those obtained by other test methods due to the effects of strain and load on the transformation. For example, transformation temperatures measured in accordance with Test Method F2082 will differ from those measured by the current standard.

5.5 The use of this test method for finished or semi-finished components without annealing (as in 8.2) shall be agreed upon between the purchaser and the supplier.

¹ This test method is under the jurisdiction of ASTM Committee F04 on Medical and Surgical Materials and Devices and is the direct responsibility of Subcommittee F04.15 on Material Test Methods.

Current edition approved Dec. 1, 2016. Published January 2017. Originally approved in 2000. Last previous edition approved in 2010 as F2004 – 05 (2010). DOI: 10.1520/F2004-16.

² For referenced ASTM standards, visit the ASTM website, www.astm.org, or contact ASTM Customer Service at service@astm.org. For Annual Book of ASTM Standards volume information, refer to the standard's Document Summary page on the ASTM website.

6. Interferences

6.1 Make sure the material to be tested is homogeneous since milligram sample quantities are used.

6.2 Take care in preparing the sample (1). Cutting and grinding can cause localized heating and/or deformation, that affect the transformation temperature. Oxidation during heat treatment can change the thermal conductance of the sample.

6.3 Set the gas flow to provide adequate thermal conductivity in the test cell.

7. Apparatus

7.1 Use a differential scanning calorimeter capable of heating and cooling at rates up to 10°C/min and of automatically recording the differential energy input between the specimen and the reference to the required sensitivity and precision.

7.2 Use sample capsules or pans composed of aluminum or other inert material of high thermal conductivity.

7.3 Use helium gas purge supply. See 10.3.1.

7.4 Use an analytical balance with a capacity of 100 mg capable of weighing to the nearest 0.1 mg.

8. Sampling

8.1 Use a sample size of 25 to 45 mg. Cut the sample to maximize surface contact with the (DSC) sample pan.

8.2 Anneal the sample at 800 to 850°C for 15 to 60 min in vacuum or inert atmosphere, or in air with adequate protection from oxidation. Rapidly cool the sample to prevent precipitation of phases that may change the transformation temperature of the alloy. For example, bare 25-45 mg samples can be air cooled after solution annealing.

8.3 Clean the sample of all foreign materials such as cutting fluid. If the sample is oxidized during heat treatment, grind, polish, or etch the sample to remove the oxide. Take care to avoid cold working the sample as this will change its thermal response. Slight oxidation is permissible but remove all heavy oxide scale.

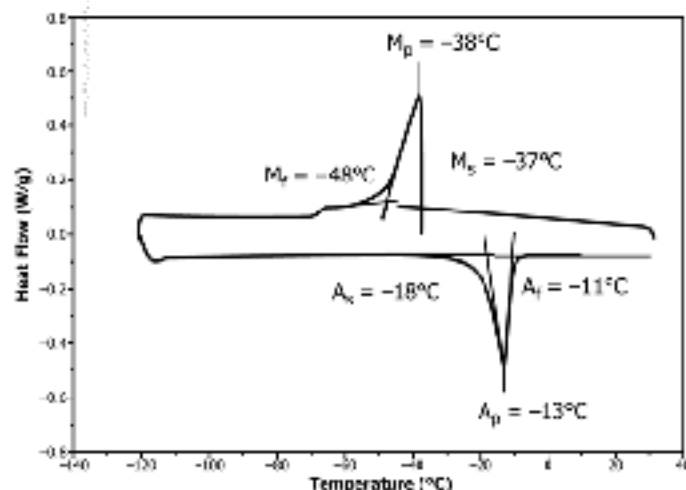


FIG. 1 DSC Curve for Nickel-Titanium (NiTi)

9. Calibration

9.1 Calibrate the temperature axis of the instrument using the same heating rate, purge gas, and flow rate as those used for analyzing the specimen in accordance with Test Method E967.

10. Procedure

10.1 Place the sample on the sample pan and place the pan on the test pedestal.

10.2 Place an empty pan on the reference pedestal.

10.3 Turn on the purge gas at a flow rate of 10 to 50 mL/min.

10.3.1 Use helium as the purge gas for the sample chamber.

10.3.2 Use a dry air, helium, or nitrogen cover gas. The dry gas shall have a dew point below the lowest temperature of the cooling cycle.

10.4 Run the cooling and heating program.

10.4.1 Use the heating and cooling rates of $10 \pm 0.5^\circ\text{C}/\text{min}$.

10.4.2 Heat the sample from room temperature to a temperature of at least $A_f + 30^\circ\text{C}$; hold at that temperature for a time sufficient to equilibrate the sample with the furnace.

10.4.3 Cool the sample to a temperature of below $M_f - 30^\circ\text{C}$; hold for a time sufficient to equilibrate the sample with the furnace. Then, heat the sample to a temperature of at least $A_f + 30^\circ\text{C}$.

10.5 *Data Acquisition*—Record the resulting curve from the heating and cooling program from $A_f + 30^\circ\text{C}$ to $M_f - 30^\circ\text{C}$.

11. Graphical Data Reduction

11.1 Draw the baselines for the cooling and heating portions of the curve as shown in Fig. 1.

11.2 Draw the tangents to the cooling and heating spikes through the inflection points as shown in Fig. 1. If a computer program is used to construct the tangents, care must be taken in locating the tangent points.

11.3 Obtain M_s , M_f , A_s , and A_f as the graphical intersection of the baseline with the extension of the line of maximum inclination of the appropriate peak of the curve as shown in Fig. 1. A_p is the peak minimum of the endothermic curve, and M_p is the peak maximum of the exothermic curve. Read A_p and M_p directly from the graph as shown in Fig. 1.

12. Report

12.1 Report the following information with the test results:

12.1.1 Complete identification and description of the material tested including the specification and lot number.

12.1.2 Description of the instrument used for the test.

12.1.3 Statement of mass, dimensions, and geometry.

12.1.4 Material for the specimen pan and temperature program.

12.1.5 Description of the temperature calibration procedure.

12.1.6 Identification of the specimen environment by gas, flow rate, purity, and composition.

12.1.7 Results of the transformation measurements using the nomenclature in accordance with Terminology F2005. Temperature results should be reported to the nearest 1°C .