

Designation: F2267 - 24

Standard Test Method for Measuring Load-Induced Subsidence of Intervertebral Body Fusion Device Under Static Axial Compression¹

This standard is issued under the fixed designation F2267; the number immediately following the designation indicates the year of original adoption or, in the case of revision, the year of last revision. A number in parentheses indicates the year of last reapproval. A superscript epsilon (ε) indicates an editorial change since the last revision or reapproval.

1. Scope

- 1.1 This test method specifies the materials and methods for the axial compressive subsidence testing of non-biologic intervertebral body fusion devices, spinal implants designed to promote arthrodesis at a given spinal motion segment.
- 1.2 This test method is intended to provide a basis for the mechanical comparison among past, present, and future non-biologic intervertebral body fusion devices. This test method is intended to enable the user to mechanically compare intervertebral body fusion devices and does not purport to provide performance standards for intervertebral body fusion devices.
- 1.3 This test method describes a static test method by specifying a load type and a specific method of applying this load. This test method is designed to allow for the comparative evaluation of intervertebral body fusion devices.
- 1.4 Guidelines are established for measuring test block deformation and determining the subsidence of intervertebral body fusion devices.
- 1.5 Since some intervertebral body fusion devices require the use of additional implants for stabilization, the testing of these types of implants may not be in accordance with the manufacturer's recommended usage.
- 1.6 *Units*—The values stated in SI units are to be regarded as the standard with the exception of angular measurements, which may be reported in terms of either degrees or radians.
- 1.7 The use of this standard may involve the operation of potentially hazardous equipment. This standard does not purport to address all of the safety concerns, if any, associated with its use. It is the responsibility of the user of this standard to establish appropriate safety, health, and environmental practices and determine the applicability of regulatory limitations prior to use.
- 1.8 This international standard was developed in accordance with internationally recognized principles on standard-

¹ This test method is under the jurisdiction of ASTM Committee F04 on Medical and Surgical Materials and Devices and is the direct responsibility of Subcommittee F04.25 on Spinal Devices.

Current edition approved April 15, 2024. Published April 2024. Originally approved in 2003. Last previous edition approved in 2022 as F2267 – 22. DOI: 10.1520/F2267-24.

ization established in the Decision on Principles for the Development of International Standards, Guides and Recommendations issued by the World Trade Organization Technical Barriers to Trade (TBT) Committee.

2. Referenced Documents

- 2.1 ASTM Standards:²
- E4 Practices for Force Calibration and Verification of Testing Machines
- E177 Practice for Use of the Terms Precision and Bias in ASTM Test Methods
- E691 Practice for Conducting an Interlaboratory Study to Determine the Precision of a Test Method
- F1839 Specification for Rigid Polyurethane Foam for Use as a Standard Material for Testing Orthopaedic Devices and Instruments
- F2077 Test Methods for Intervertebral Body Fusion Devices

3. Terminology

- 3.1 All subsidence testing terminology is consistent with the referenced standards above, unless otherwise stated.
 - 3.2 Definitions:
- 3.2.1 coordinate system/axes—three orthogonal axes are defined following a right-handed Cartesian coordinate system (Fig. 4). The XY plane bisects the sagittal plane between the superior and inferior surfaces that are intended to simulate the adjacent vertebral end plates. The positive Z-axis is to be directed superiorly. Force components parallel to the XY plane are shear components of loading. The compressive axial force is defined to be the component in the negative Z direction. Torsional load is defined to be the component of moment about the Z-axis.
- 3.2.1.1 *origin*—the center of the coordinate system is located at the center of rotation of the testing fixture.
- 3.2.1.2 *X-axis*—the positive X-axis is a global fixed axis relative to the testing machine's stationary base and is to be directed anteriorly relative to the specimen's initial unloaded position.

² For referenced ASTM standards, visit the ASTM website, www.astm.org, or contact ASTM Customer Service at service@astm.org. For *Annual Book of ASTM Standards* volume information, refer to the standard's Document Summary page on the ASTM website.

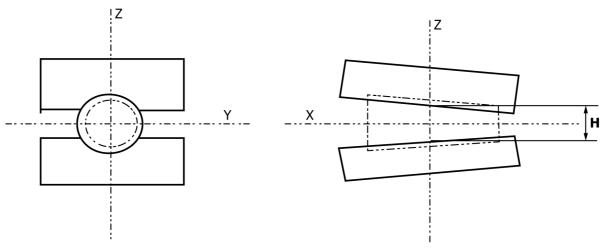


FIG. 1 Intradiscal Height Diagram

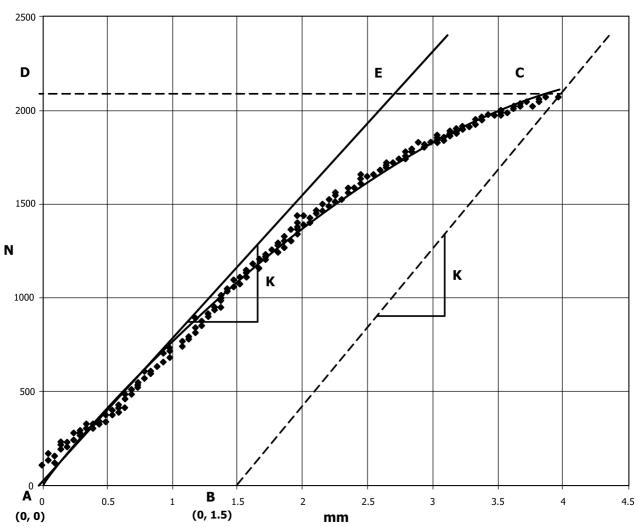


FIG. 2 Typical Load-Displacement Curve with 1.5 mm (Thoracic Device) Offset for Polyurethane Foam Test Blocks

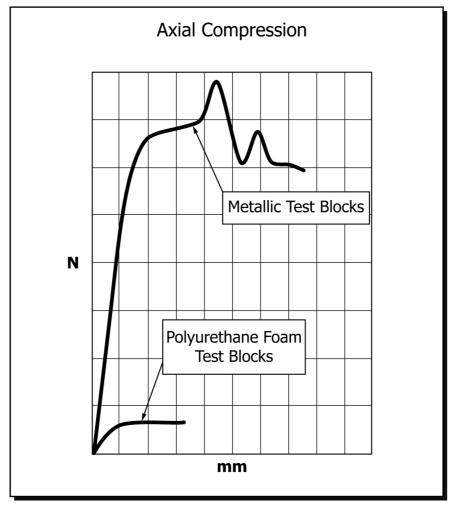


FIG. 3 Typical Load-Displacement Plot Comparison for Test Specimens in Metallic and Polyurethane Test Blocks

- 3.2.1.3 *Y-axis*—the positive Y-axis is a global fixed axis relative to the testing machine's stationary base and is directed laterally relative to the specimen's initial unloaded position.
- 3.2.1.4 *Z-axis*—the positive *Z*-axis is a global fixed axis relative to the testing machine's stationary base and is to be directed superiorly relative to the specimen's initial unloaded position.
- 3.2.2 *ideal insertion location*—the implant location with respect to the simulated inferior and superior vertebral bodies (polyurethane) dictated by the type, design, and manufacturer's surgical installation instructions.
- 3.2.3 intended method of application—intervertebral body fusion device assemblies may contain different types of stabilizing anchors such as threads, spikes, and knurled surfaces. Each type of anchor has an intended method of application or attachment to the spine.
- 3.2.4 intended spinal location—the anatomic region of the spine intended for the intervertebral body fusion device assembly. Intervertebral body fusion device assemblies may be designed and developed for specific regions of the spine such

- as the lumbar, thoracic, and cervical spine. Also, there exist different anatomical potential surgical approaches, which will result in different implant orientation at different levels of the spine.
- 3.2.5 *intervertebral subsidence*—the process of a vertebral body cavitating or sinking around an implanted intervertebral body fusion device resulting in the loss of intradiscal height.
- 3.2.6 *intradiscal height*—the straight-line distance along the Z-axis between the unaltered simulated vertebral bodies. See Fig. 1.
- 3.2.7 *load point*—the point through which the resultant force on the intervertebral device passes; that is, the geometric center of the superior fixture's sphere (see Fig. 4).
- 3.2.8 offset displacement—offset on the displacement axis equal to 1 mm for cervical disc devices, 1.5 mm for thoracic devices, and 2 mm for lumbar devices (see distance AB in Fig. 2)
- 3.2.9 *simulated vertebral bodies*—the component of the test apparatus for mounting the intervertebral body fusion device.